
TRICHOTRACK LITE

Installation Manual

Version 1.0

Version	Date	Author	Change Description
1.0	26 Feb 2024	Xiaoxing Liu	Document created



1 Installation Requirements

Before installing TrichoTrack Lite, ensure that your device meets the following requirements:

1.1 System Requirements

- Tablet: iPad 10th generation or above, iPad Air 5th generation or above.
- Operation System: iPadOS 17.3 or above.
- Web Browser: Safari of the iPadOS.

1.2 Internet Connectivity

TrichoTrack Lite requires an active internet connection for updates and data synchronization. Make sure your device has reliable internet access.

1.3 Digital Microscope Compatibility

TrichoTrack Lite requires a digital microscope. Currently, it only supports the below model:


- I. Plugable 3691 60x-250x Digital Microscope
 - With hybrid 2-in-1 USB 2.0 and USB-C cable.

Please ensure you have a microscope of the above model with the USB-C connection for accurate scalp scan.


2 Installation Instructions

Start to install TrichoTrack Lite by opening the Safari web browser.


1. Go to the website <https://thriving-lokum-87d5f6.netlify.app/>. Click the square icon on the right of the address bar to initiate the App installation.

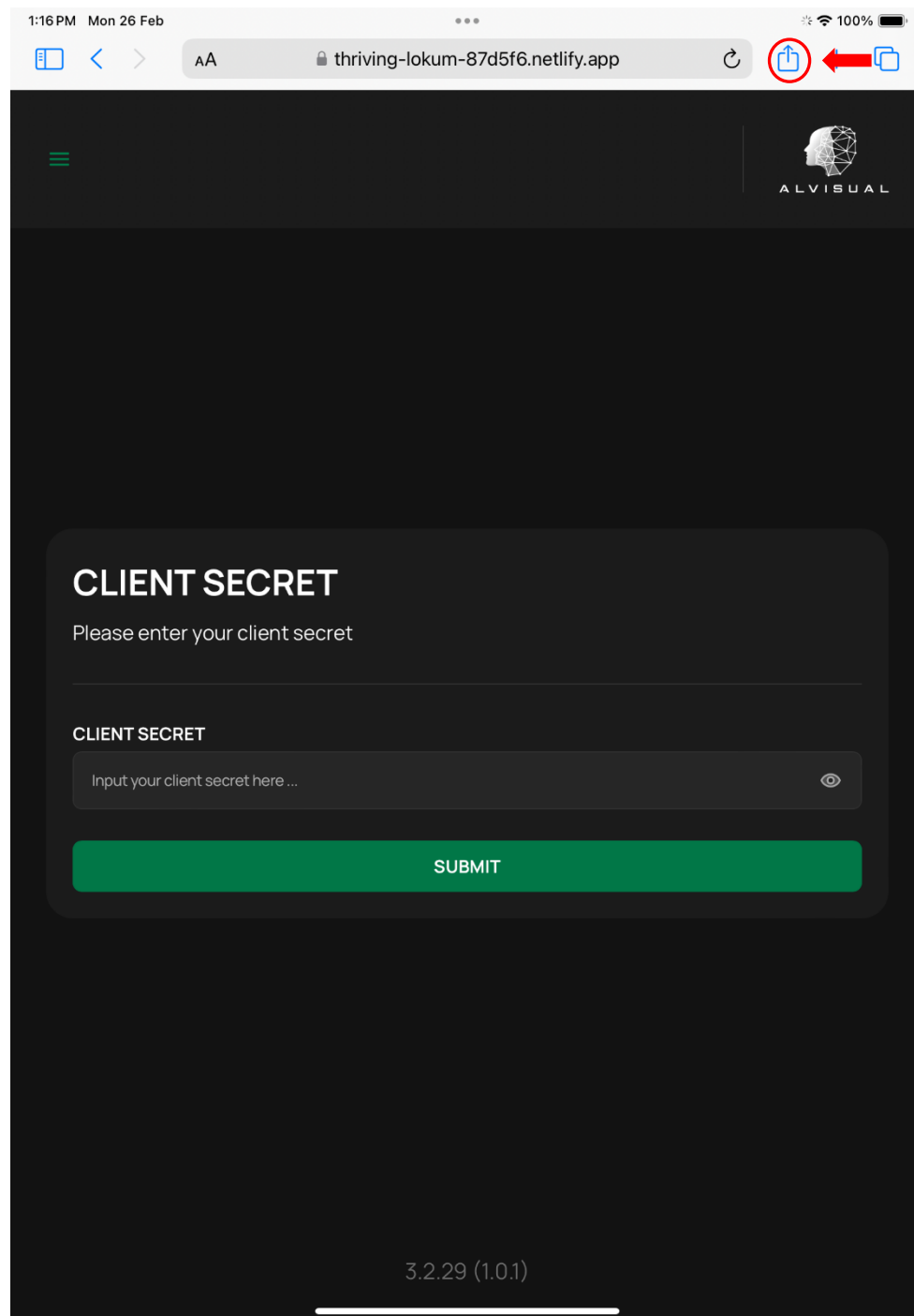

2. In the drop-down menu, click 'Add to Home Screen'.

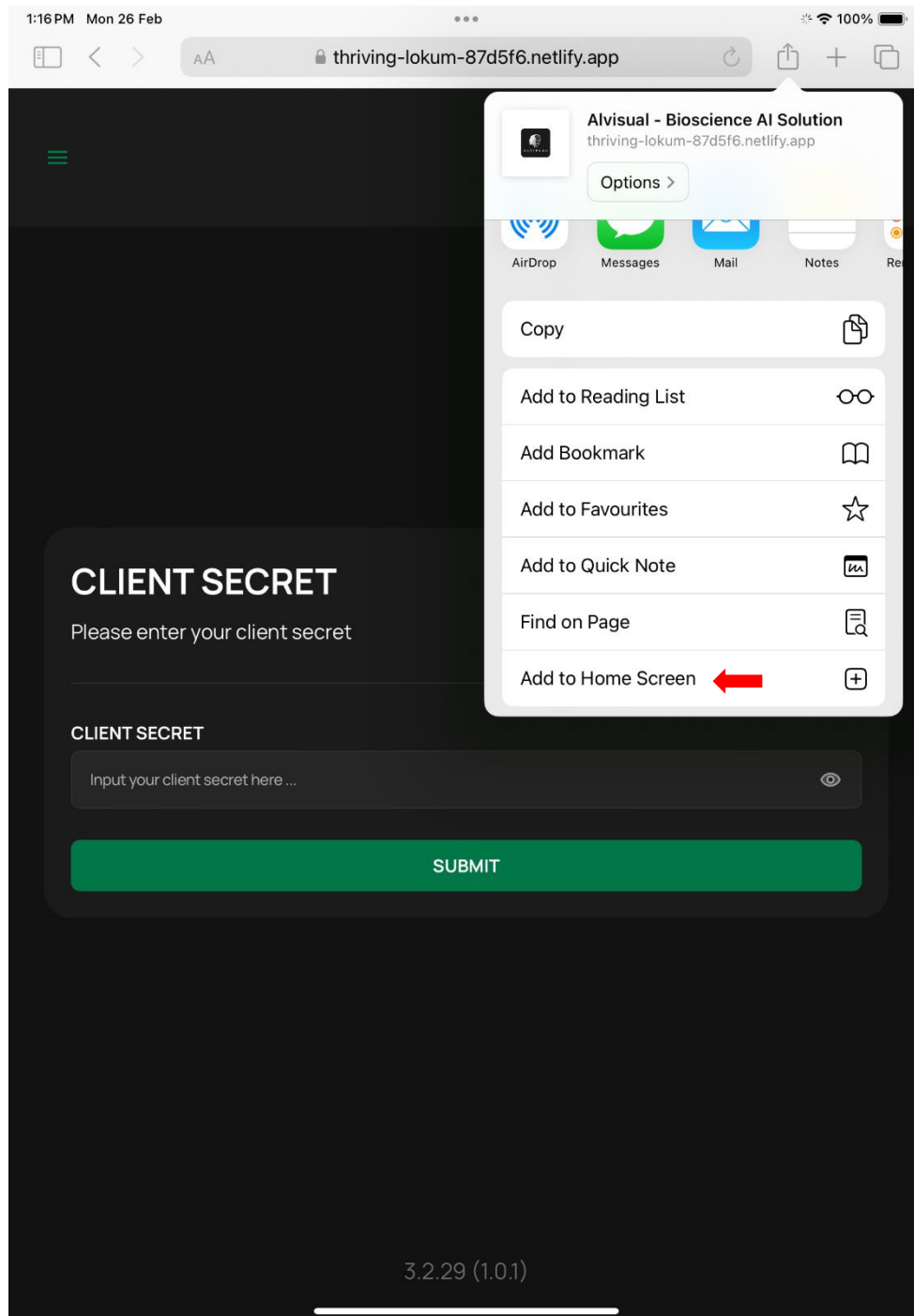

3. A dialog window will appear to allow you to add the App to the iPad home screen. Click 'Add' button.

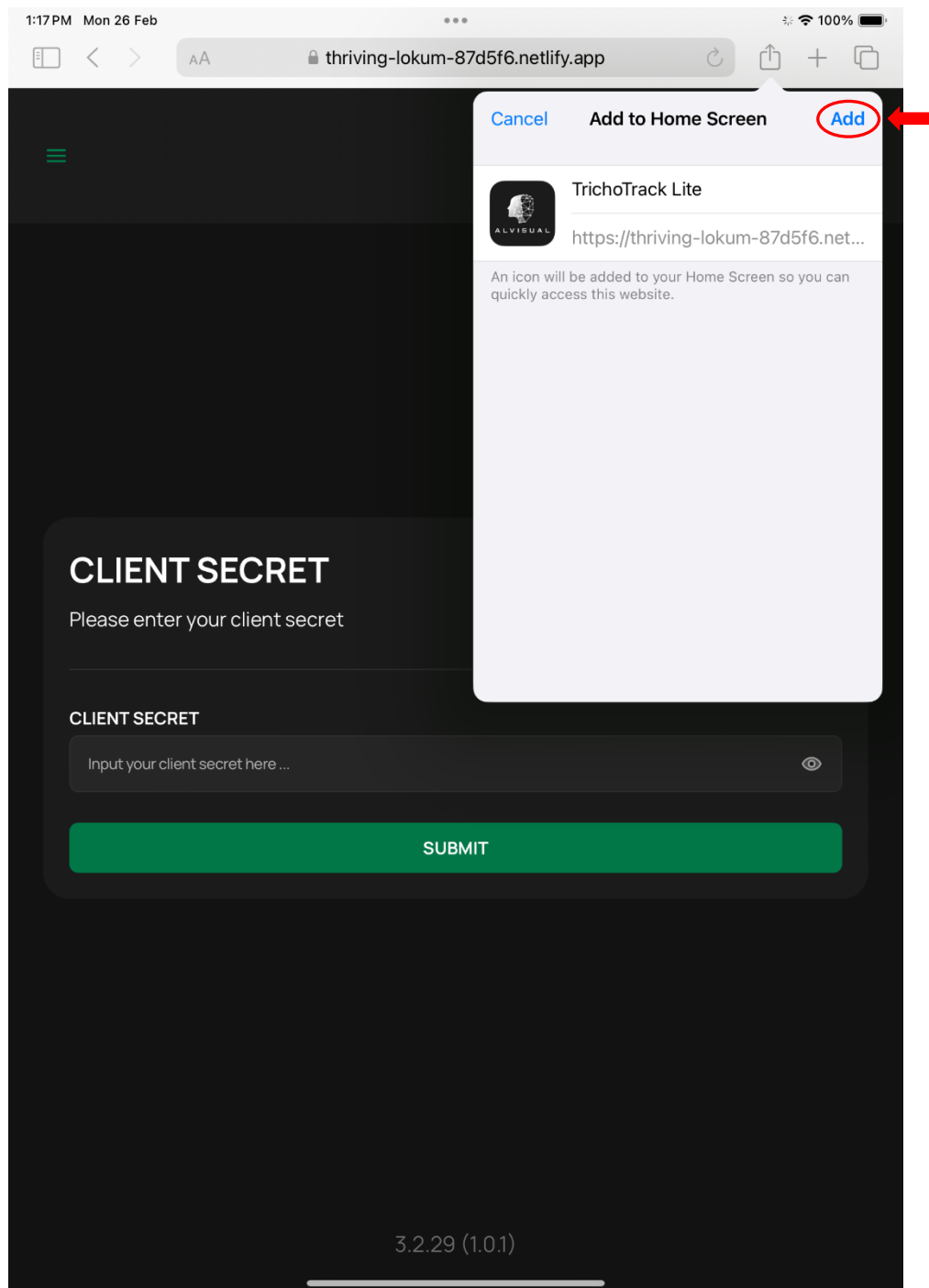



The installation process is now completed. You'll find a new app icon of 'TrichoTrack Lite' on your iPad home screen. To start TrichoTrack Lite, simply double-click the icon. For the best experience, run it in **portrait mode**.

